
Panasonic



Operating Instructions
(Household) Oral Irrigator


Model No. EW1611

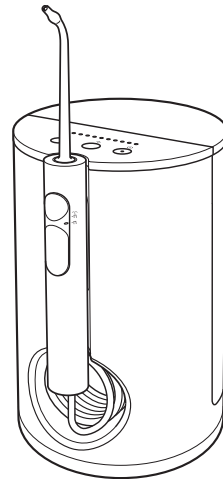CE

Ⓒ English	5	Ⓒ Nederlands	53	Ⓒ Česky	101
Ⓒ Deutsch	17	Ⓒ Español	65	Ⓒ Slovensky	113
Ⓒ Français	29	Ⓒ Português	77	Ⓒ Magyar	125
Ⓒ Italiano	41	Ⓒ Polski	89	Ⓒ Română	137

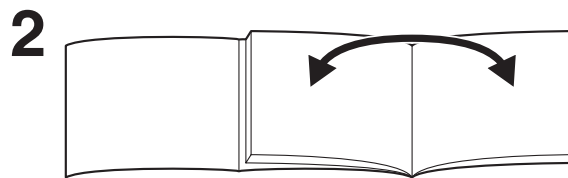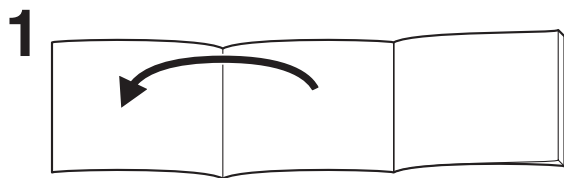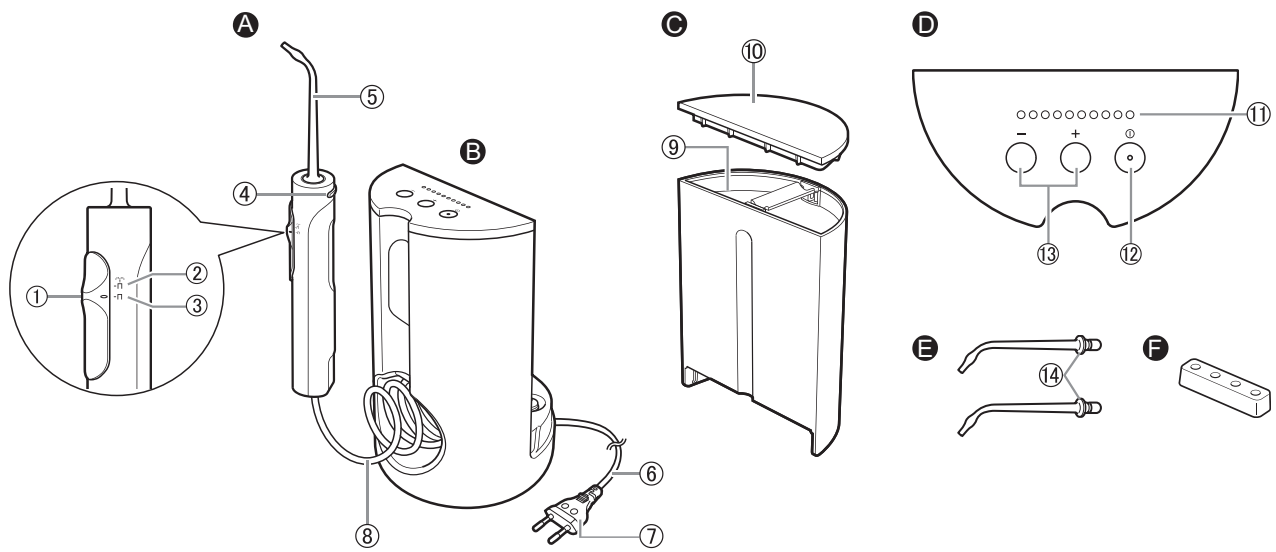Panasonic

Operating Instructions
(Household) Oral Irrigator

Model No. EW1611

Contents

Safety precautions	8	Troubleshooting	14
Intended use	10	Specifications	15
Parts identification	11	Changing the nozzle	
How to use	11	(Sold separately).....	15
After use	13		

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

5

Warning

- Do not immerse the appliance in water or wash it with water.
Doing so may cause electric shock or fire due to a short circuit.
- Always unplug the power plug from a household outlet when cleaning the power plug.
Failure to do so may cause electric shock or injury.
- The supply cord cannot be replaced. If the cord is damaged, the appliance should be scrapped.
Failure to do so may cause an accident or injury.
- Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.

- Oral irrigators shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions.
- Children shall not play with the appliance.

Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

 WARNING	Denotes a potential hazard that could result in serious injury or death.
 CAUTION	Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe the type of instructions to be observed.

	This symbol is used to alert users to a specific operating procedure that must not be performed.
	This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

WARNING

► Preventing accidents

Do not store within the reach of children or infants.

Do not let them use it.

- Doing so may cause an accident or injury due to accidental ingestion of removable parts.

Properly dispose of the appliance package, and keep it away from children or infants.

- Failure to do so may cause suffocation or other accidents to children or infants.

► Power supply

Do not connect or disconnect the power plug to a household outlet with a wet hand.

- Doing so may cause electric shock or injury.

Do not use the appliance near a bath tub or inside a shower room.

- Doing so may cause electric shock or fire.

Do not use when the cord or the power plug is damaged or when the fitting into the household outlet is loose.

Do not damage, modify, forcefully bend, pull, twist, or bundle the cord.

Also, do not place anything heavy on or pinch the cord.

- Doing so may cause electric shock or fire due to a short circuit.

WARNING

Do not use in a way exceeding the rating of the household outlet or the wiring.

- Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.

Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the appliance.

Fully insert the power plug.

- Failure to do so may cause fire, electric shock, or injury.

► **In case of an abnormality or malfunction**

Immediately stop using and remove the power plug if there is an abnormality or malfunction.

- Failure to do so may cause fire, electric shock, or injury.

<Abnormality or malfunction cases>

- **The main unit or cord is deformed or abnormally hot.**
 - **The main unit or cord smells of burning.**
 - **There is an abnormal sound during use of the main unit or cord.**
- Immediately request inspection or repair at an authorized service centre.

► **This product**

Do not drop or subject to shock.

- Doing so may cause injury.

Do not modify, disassemble or repair.

- Doing so may cause injury due to fire or malfunction.
- Contact an authorized service centre for repair.

► **Cleaning**

Regularly clean the power plug to prevent dust from accumulating (once every half a year).

- Failure to do so may cause fire due to insulation failure caused by humidity.

Disconnect the power plug and wipe with a dry cloth.

CAUTION

► **Preventing damage to teeth and gums**

People who are unable to operate should not use this appliance.

People with weak oral senses should not use this appliance.

- Doing so may cause injury or damage to teeth and gums.

People who may have periodontal disease, have teeth being treated, or concerned about symptoms within their mouth should consult a dentist before use.

- Failure to do so may cause injury or damage to teeth and gums.

GB

D

F

I

NL

E

P

PL

CZ

SK

H

RO

CAUTION

► This product

Never use for a purpose other than oral cleaning.

-  - Doing so may cause an accident or injury when the stream of water is pointed toward the eyes, nose, ears, or throat.

-  **Do not fill the tank with water warmer than 40 °C.**
- Doing so may cause burns.

-  **Do not share nozzle with your family or other people.**
- Doing so may cause infection or inflammation.

► Note the following precautions

Do not allow metal objects or trash to adhere to the power plug.

-  - Doing so may cause electric shock or fire due to a short circuit.

Do not wrap the cord around the appliance when storing.

-  - Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.

Disconnect the power plug from the household outlet when not using for extended period of time.

-  - Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

Disconnect the power plug by holding onto it instead of the cord.

-  - Failure to do so may cause electric shock or fire due to a short circuit.

After use, always care for and clean the appliance. (See page 13.)

-  - Failure to do so may cause damage to your health if used without keeping it sanitized due to propagation of mold and bacteria.

Intended use

► Notice

- Do not put salt, mouthwash or other chemicals into the water tank. Doing so may cause malfunction.
- Do not dry the main unit with a dryer or a fan heater. Doing so may cause malfunction or deformation of the parts.
- Do not wipe with thinners, benzine, alcohol etc. Doing so may cause malfunction, cracking or discoloration of the parts.
- Store the appliance in a place with low humidity after use. Leaving the appliance in a bathroom may cause malfunction.

Parts identification

A Handle

- ① Stop/Jet switch
- ② Jet mode
- ③ Stop
- ④ Nozzle release button
- ⑤ Nozzle

B Main body

- ⑥ Cord
- ⑦ Power plug
- ⑧ Suction hose

C Water tank

- ⑨ Upper water level limit
- ⑩ Water tank cap

D Operation panel

- ⑪ Water pressure level LED indicator (1~10 levels)

- ⑫ Power switch
- ⑬ Water pressure regulator switch

E Nozzle (x2)

- The nozzle is a consumable.

- ⑭ Identification ring (white & blue)

Accessory


F Nozzle stand

There may be moisture inside the main unit (water tank and suction hose) left from the distilled water used for product testing.

How to use

On sink or shelf

- ① Place the main unit on a flat and stable surface.
- ② Plug in the power plug.
 - Convenient dual power cord outlet able to pull out the cord, both right and left.

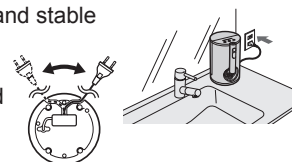1. Open the water tank cap, remove the water tank, and fill it with water.


- Always use new tap water or lukewarm water.
- If the water feels cold, use lukewarm water with a temperature of 40 °C or less.
- When filling water, be sure not to exceed the upper water level limit.

2. Set the tank to the main unit firmly, and attach the water tank cap.

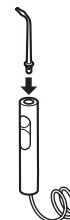
3. Insert the nozzle firmly into the handle. Make sure fixed firmly.



11

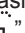4. Select the water pressure level.

- ① Press the power switch.
(The water pressure level LED indicator lights on.)
 - The product will automatically stop operation after about 10 minutes.
- ② With the nozzle facing the basin, set the Stop/Jet switch to “” (Jet mode) to jet water.
- ③ Press “+” or “-” on the water pressure regulator switch until the desired flow intensity is reached.
 - Keep pressing “+” or “-” will continuously change the water pressure.
- ④ Set the Stop/Jet switch to “” (Stop) to stop jetting temporarily.
 - Do not leave the product in the “” (Stop) state for more than 1 minute. (May cause malfunction.)


Stop/Jet switch

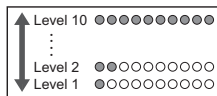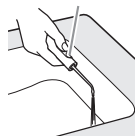5. Place the nozzle in your mouth, close your mouth gently, and set the Stop/Jet switch to “” (Jet mode).



If you want to pause jetting, set the Stop/Jet switch to “” (Stop).

Raise your elbow — Keep head over the basin

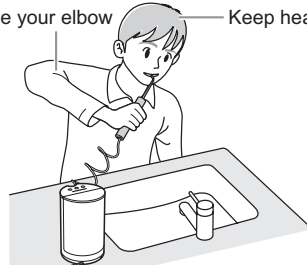6. After use, set the Stop/Jet switch to “” (Stop) and turn off the power switch.

- The water pressure level LED indicator flashes. (The indicator will be off after 1 minute.)
- The last selected water pressure level will be remembered, from which the oral irrigator will start next time.

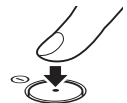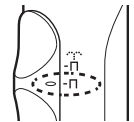► **Notes**

- Do not push the tip of the nozzle too hard against the teeth or gums.
- Since plaque adhering to the surface of teeth cannot be removed only with water jet, brushing is always necessary. It is recommended to clean inside your mouth with water jet after brushing your teeth.
- Do not operate the main unit when empty other than when pouring out any remaining water after use. (Operating the main unit without water may result in a malfunction.)

Cleaning the teeth

Interdental cleaning

- Direct the water jet onto the space between the teeth.
- Direct from the front and rear sides of the teeth.

Periodontal pocket cleaning

- Aim the water jet at the space between the teeth and gum.
- Slowly shift the water jet along the gum.
- Clean the front and back of the teeth and behind the molars as well.

Gum care

- Stimulate the gum.
- If you are concerned about your gum, tenderly stimulate it with gentle flow of water.


After use

1. Take off the water tank cap and throw away the remaining water and wash.

- Change the water every time you use the irrigator.
- The water tank can be washed in a dish washer.
- The water tank cap cannot be cleaned in a dish washer.

2. Hold the nozzle over the basin, before the Stop/Jet switch to “” (Jet mode).

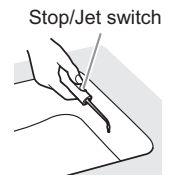
3. Turn on the power switch to drain the remaining water from the main unit.

- Run the irrigator all the way until there is no water coming out, and turn off the power switch.


4. Disconnect the power plug from the outlet.

5. Push the release button on the body to release the nozzle, then wash it. Keep the nozzle.

- The nozzle stand cannot be cleaned in a dish washer.

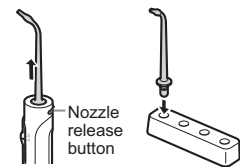6. Clean the handle, and put it back into the storage slot.

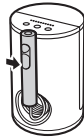When cleaning the nozzle and handle

- Do not use lukewarm water exceeding 50 °C.
- Use neutral kitchen detergent and rinse well with water to ensure no detergent remains after cleaning.
- Wipe away the water with a towel or the like, keep the product in the dry place.
- The nozzle cannot be cleaned in a dish washer.

When wiping stains away on the main unit with a cloth

- Never use thinners, benzine, or alcohol. (Doing so may cause malfunction, cracking or discoloration of the parts.)
- Do not wash the main body with water.

Troubleshooting

Problem	Possible cause	Action
No water comes out, or water jet power down.	<ul style="list-style-type: none"> ▶ The nozzle is clogged. ▶ The outlet of nozzle is squashed. 	▶ Replace the nozzle.
The nozzle comes off during use.	▶ The nozzle is deformed.	
Motor sound changes after the Stop/Jet switch is set to "Π" (Stop).	▶ This is not an abnormality.	
Water jet stops during using.	▶ Takes 10 minutes after switch on.	▶ This is not an abnormality. (See page 12.)
	▶ The self-protection function is activated due to a thunderstorm etc.	▶ Turn on the power switch again.
The water pressure level LED indicator still flashes after the power switch is turned off.	▶ This is not an abnormality. The indicator will be off after 1 minute of turning off.	

If the problems still cannot be solved, contact the store where you purchased the unit or a service centre authorized by Panasonic for repair.

Specifications

Power Source	See the name plate on the product (Automatic voltage conversion)
Power consumption	18 W
Water tank volume	Approx. 600 mL (the water pressure level "10", 115 seconds)
Dimensions	Main unit: 280 (H) × 125 (W) × 140 (D) mm (Including nozzle) 200 (H) × 125 (W) × 140 (D) mm (Not including nozzle)
Suction hose length	Approx. 1100 mm
Mass	Main unit: Approx. 940 g (Including nozzle)
Airborne Acoustical Noise	82 (dB (A) re 1 pW)

This product is intended for household use only.

Changing the nozzle (Sold separately)

- If the tip of nozzle gets deformed, please replace with a new nozzle.
- The nozzle is a consumable. The nozzle should be replaced with a new one every 6 months or so due to hygiene reasons, even if its shape is unchanged.

New nozzle (tip)

Deformed (worn) nozzle (tip)

Part number	Replacement Nozzle for EW1611	
WEW0982	Nozzle	

Disposal of Old Equipment
Only for European Union and countries with recycling systems

This symbol on the products, packaging, and/or accompanying documents means that used electrical and electronic products must not be mixed with general household waste. For proper treatment, recovery and recycling of old products, please take them to applicable collection points in accordance with your national legislation.

By disposing of them correctly, you will help to save valuable resources and prevent any potential negative effects on human health and the environment.

For more information about collection and recycling, please contact your local municipality.

Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

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© Panasonic Corporation 2017
Panasonic Corporation